# smovey neuroVIBE

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...to promote and activate the functions of the  
musculoskeletal system and the nervous system

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Regular exercise and physical activity have a positive impact on your state of health and can even replace and/or reduce the medication you take. The benefits of exercise have been studied intensively for more than 30 years and the research has shown that in addition to smoking a lack of exercise and sitting for too long are the major risk factors for various illnesses (cf. Löllgen et al, 2018). People with a healthy lifestyle contribute around 50% to their health with regular physical activity (cf. Löllgen et al, 2018).

It is amazing that despite these findings, exercise has not gained much more in importance; if it were given the status it deserves, a whole number of illnesses could be avoided by way of prevention. Plus, it would also be possible to help treat illnesses such as high blood pressure, osteoporosis, cancers, neurological disorders, musculoskeletal disorders etc. with exercise (vgl. Löllgen et al, 2018).

This is where smovey's philosophy comes into its own: regular exercise that is fun based on the principle of swing, move and smile. Because only if physical activity is enjoyable are we human beings capable of practising it in the long term and thus take responsibility for improving our everyday life and health.

With the smovey neuroVIBE we have created an exercise device that can be used in a variety of ways for both endurance and strength training. Plus, it uses visual, acoustic and haptic stimuli along with existing ones such as vibration in a unique way to support the functions of the musculoskeletal system and the activity of your nervous system.

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Eure Marion Reisinger und Team

The smovey neuroVIBE is an excellent option for people who wish to include some exercise and physical activity in their lives again after a period of inactivity or who are adversely affected by physical limitations. It can be used at any level of fitness and the effects of even short fitness sessions can be felt very rapidly thanks to the intensive stimuli.

Physical activity as medicine (Löllgen et al, 2018) is coming a little closer with the smovey neuroVIBE and should help us to prescribe exercise preferably for ourselves every day - for a life of health and fitness, ideally by way of prevention.

Have lots of fun doing regular exercise with the smovey neuroVIBE.

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Source:

Prof Herbert Löllgen, Remscheid, Dr Jürgen Wismach, Berlin and Prof Norbert Bachl, Vienna Physical activity as medicine, Drug prescription in practice, volume 45, issue 3, July 2018

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**smovey neuroVIBE** is an oscillating ring system consisting of a spiral tube, 3 or 4 steel balls (depending on the model) soft damping elements and 2 handle systems. Oscillation of the rings forwards, backwards or in a circular formation prompts the steel balls to move as free mass, generating the vibration that is typical of smovey by way of their rolling motion across the grooves and ridges of the spiral. This vibration, which occurs within a frequency range that is in consonance with the human organism, is immediately palpable in the palms of the hands.

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
### **smovey neuroVIBE - small/light**

A smaller and lighter model with 3 balls and a soft foam rubber grip. Excellent as a device for children and when used by individuals with small hands, for people who start again after a long break in training or for people with compromised strength.

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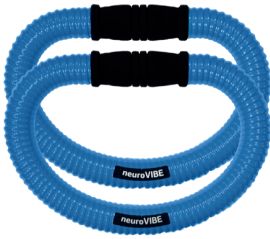### **smovey neuroVIBE - Classic**

The classic smovey neuroVIBE, featuring 4 balls and a soft foam rubber grip, is ideal when used as a general preventive device by adults and older persons in the field of movement and the targeted development and improvement of the functions of the musculoskeletal and nervous system.

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### **smovey neuroVIBE - Aqua**

The waterproof design of smovey neuroVIBE - Aqua also makes it suitable for use in all types of water settings. It features an ergonomic, firmer and quick-drying grip and 4 balls. Very easy to disinfect with standard disinfectants.

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All **smovey neuroVIBE-models** come equipped with the soft damping system, which means that the steel balls turn smoothly within the spiral tube, thus ensuring a very gentle and yet effective training session irrespective of the age of the user!

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*Stimulate the hand reflex points according to the model of Traditional Chinese Medicine (TCM).*

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*Various frequency ranges result from the size of the components and have a positive effect on the body. (TCM).*

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*The actual weight is 1.1lb at standstill. The centrifugal forces caused by the dynamic motion from swinging can increase this weight many times over.*

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*Balance exercises refine the interaction of your muscles which tends to be beneficial as you deal with sudden movements that could cause a fall.*

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*If used in certain cross motion patterns your smovey neuroVIBE activates both, the right and the left sides of your brain.*

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*Whether rain, snow or shin, at the gym or at home, due to its multifunctionality your smovey neuroVIBE is an excellent training and fitness tool.*

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*Everybody loves the outdoors. Take your smovey neuroVIBE with you to the great outdoors, the beach and the mountains and enjoy movement and active recreation.*

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
*Improve your endurance, strength and flexibility. Relieve pressure on your spine, joints, ligaments, tendons, and have great fun!*

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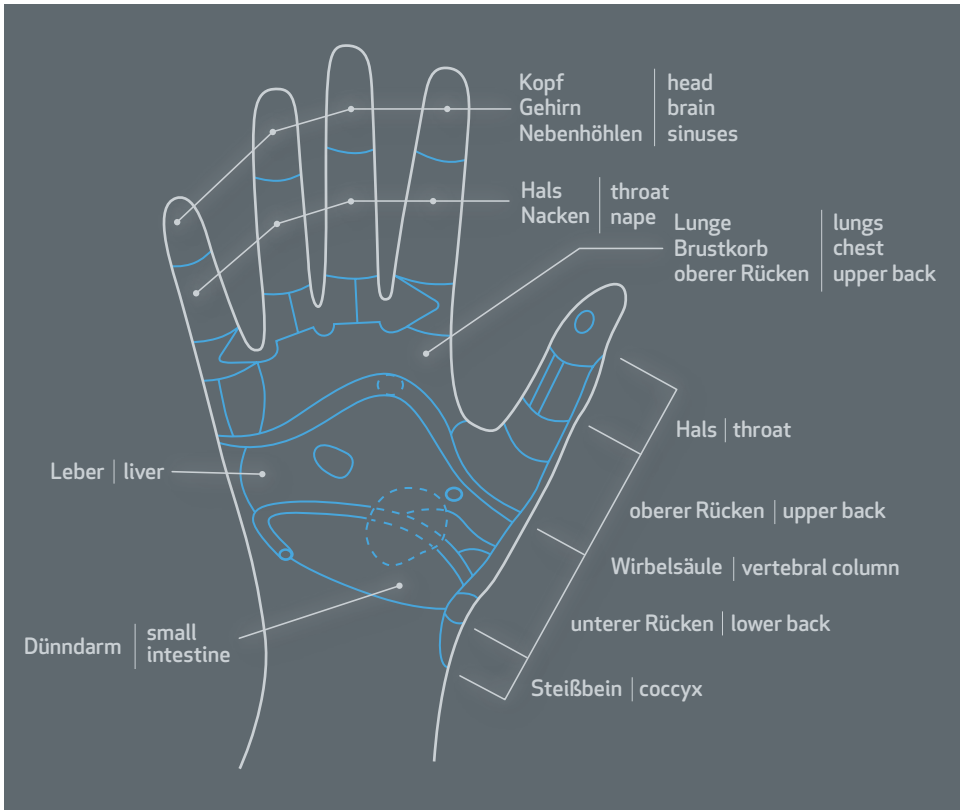The smovey neuroVIBE are produced in the town of Weyer in Oberösterreich in a „social workplace“ by hand by people with a light physical or mental limitation. After every production step a quality control check is made, and that is why smovey is a unique product: „smovey - the Original“ and „smovey - handmade in Austria“.

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Johann “SalzHans” Salzwimmer is the inventor of smovey, the vibroswing system. He stood firm in life, when he got the diagnosis “Parkinson” in the mid of the 1990s. He accepted his illness, but began, familiar with the fast development of the disease, to develop a method that would help him in stabilizing his state of health.

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smovey neuroVIBE are characterized by the particularly intensive vibration and the soft damping elements, which make the rings particularly suitable for short, intensive units of movement. Disabled people also have a reliable device here, which activates the nervous system in a variety of ways and promotes the muscles and the musculoskeletal system. Due to the different sizes, the requirements and needs of the users can be explicitly addressed and there is exactly the right model available for everyone: for support in the area of health prevention through movement for both children and adults, for training in water, for special purposes Promotion and activation of the functions of the nervous system and musculoskeletal system in people with special needs, regardless of age and level of training!

The soft damping system gently cushions the balls and enables people with disabilities to exercise gently, comfortably and effectively.

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Exercise promotes physical and mental development even at a young age. „Exercise is the gateway to learning“, as the saying goes. Children generally have a natural urge to indulge in exercise; this can be encouraged by way of the smaller smovey neuroVIBE small/light which is great fun to use. Whether combined with music, obstacle courses or dexterity exercises, whether outdoors or in the classroom - in a group setting, exercising with smovey neuroVIBE quite simply puts participants in a good mood and, in the long term, promotes powers of concentration, balance, perception and coordination. The rolling action of smovey neuroVIBE and its relaxing effect is also appreciated by children.

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Exercise in general, but particularly when facilitated by smovey neuroVIBE has a special role to play in the preventive setting. A long, independently led life requires that core competencies, such as muscular strength, balance (fall prophylaxis), coordination, intensive walking, relaxation, etc., are regularly trained in order to preserve them.

These areas are presented and the basic exercises explained on the enclosed smovey neuroVIBE training chart. The exercises can be done sitting and / or standing. Further exercises can be accessed via a link or QR code which you find at the trainingschart in the box

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**IMPORTANT:** smovey neuroVIBE is an exercise device that can positively impact the state of your health. Before a training session with smovey neuroVIBE, we strongly recommend that you are examined by a trusted physician. If you encounter problems or experience pain when using smovey neuroVIBE, a physician should be consulted at the earliest and the situation medically evaluated. Exercise and training sessions incorporating smovey neuroVIBE are conducted exclusively at your own risk. smovey GmbH assumes no liability whatsoever for any injuries or damage that might occur.

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The smoveyINVENTOR Johann 'Salzhans' Salzwimmer has suffered from Parkinson's disease for over 25 years, and he exercises using his rings 3-5 times every day, each session lasting around 5 minutes. These are some of the exercises he favours:

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The movement with smovey neuroVIBE is also particularly suitable for people who cannot stand or who have too little stability for the exercise units when standing. When you roll your legs, the vibrations can sometimes be felt in the feet and activate them. The upper body can also be trained and activated by swinging and moving the smovey neuroVIBE.

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With smovey neuroVIBE, it is possible to activate people of all ages and with different health conditions for a joint training. This promotes not only the social aspect but makes the training with smovey neuroVIBE absolutely unique.

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
**Important:**  
smovey neuroVIBE is a movement device that, with regular use, can activate and improve the functions of the nervous system and the musculoskeletal system. In the event of uncertainty or to clarify the individual situation, we urgently recommend coordinating the use of the smovey neuroVIBE rings with the treating doctor and / or therapist you trust! The smovey neuroVIBE rings are used at your own risk, smovey GmbH assumes no liability for any complaints, injuries or damage.

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Please read the following safety information carefully before you start using the smovey neuroVIBE! Follow these guidelines when using the smovey neuroVIBE - rings to ensure safety. If the smoveys are not treated or managed properly, it could lead to unforeseen malfunctions caused by fault of the user. The smovey GmbH assumes no liability for injuries, or damage!

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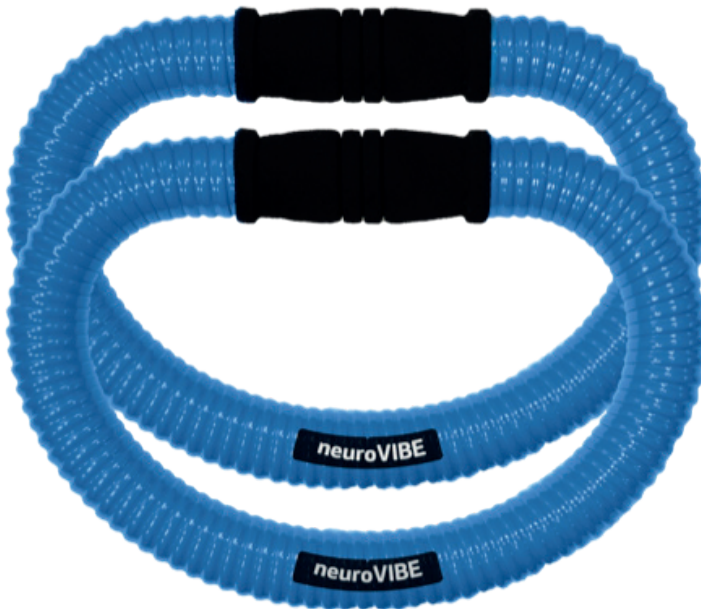
Do not expose the smovey neuroVIBE to extremely high temperatures. Avoid direct sunlight, for example, in the car or on window sills.

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Do not expose the smovey neuroVIBE to extremely low temperatures

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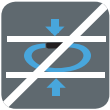Don't bend the smovey-ring! This avoids kinks in the hose. These kinks will have a direct impact on the gradient-properties of the steel balls.

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smovey neuroVIBE rings are not food and not suitable for consumption.

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Never touch or carve the smovey neuroVIBE with a sharp or pointed object (for example: knife, screw-driver, ...)! This can cause irreparable mechanical damages.

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Keep the smovey neuroVIBE rings away from open flames or burning tobacco products.

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Avoid contact and cleaning of the smovey neuroVIBE with inappropriate, corrosive or toxic liquids of gases.

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Not suitable for children under the age of 3 years.

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It is recommended, when using smovey neuroVIBE rings to wear suitable safety glasses.

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smovey GmbH | Weinzierlstraße 2a | A-4400 Steyr | +43 7252 / 716 10